



GLASS EXPANSION
Quality By Design

Instructions for

Nexus Universal Nebulizer Connection Kit

P/N FT-16-8-X



General Safety

- It is recommended you read and understand this procedure before starting.
 - It is advised you wear safety gloves and glasses when handling glass.
-


Disclaimer: Glass Expansion supplies products to customers with the assumption that only qualified individuals who are properly trained in the appropriate procedures will use our products. It is our customer's entire responsibility to safely store, handle and use the products ordered from Glass Expansion. We assume the requisitioner to be experienced, competent, safety conscious professionals. Glass Expansion assumes no liability for any damage or injuries -physical, emotional or mental, arising out of persons incorrectly or inappropriately using or handling our products. We also assume no liability for any damage to an instrument or product from incorrect installation, handling or use of the product by the user.

Instructions


1. Using a sharp razor or scalpel, cut the 1/16" sample tube to desired length - The cut end must be square and burr-free.

Caution! Be careful when handling and storing sharp objects.

2. Slide the connector over the cut end of the tube.
3. Slide 1/8" cap nut over the nebulizer sample-inlet stem. Ensure that it is positioned as shown below.
4. Slip the DualSeal ferrule over the nebulizer sample-inlet stem until it comes to a positive stop.

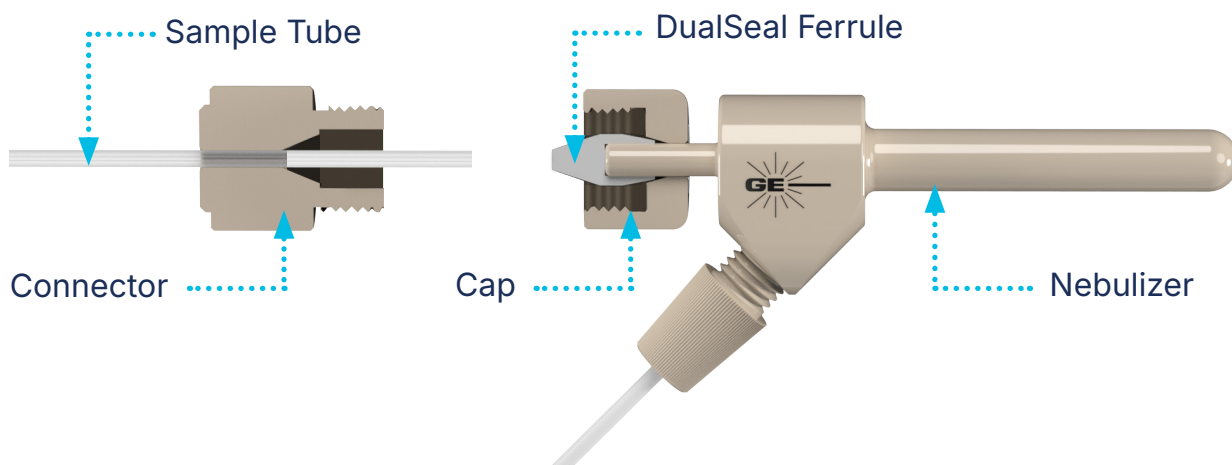

5. Holding the nebulizer and cap in one hand, push the sample tube into the DualSeal ferrule until it stops, then thread the cap onto the connector and tighten finger-tight.
6. When finished, gently pull on the sample tube to ensure a firm connection has been made.

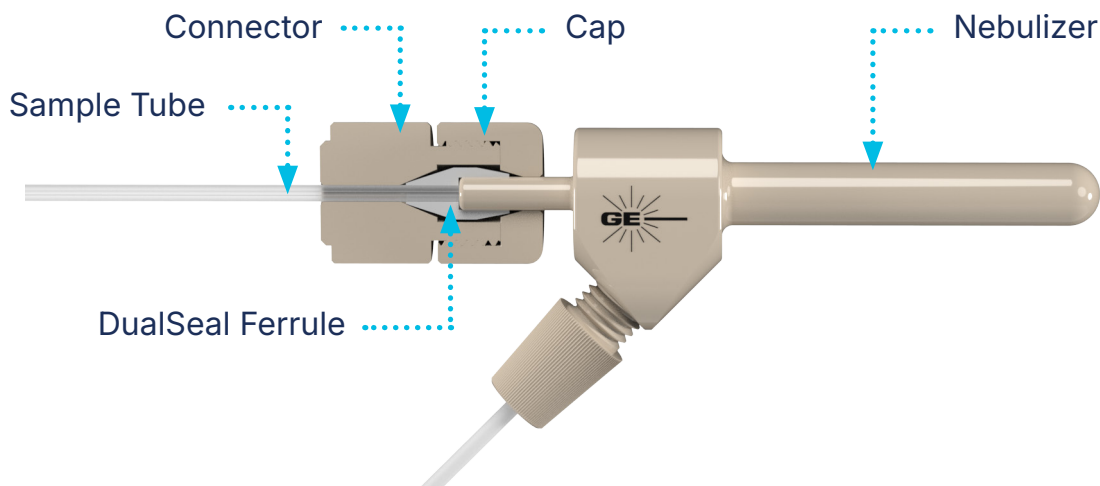

Caution! It is advised you wear safety gloves and glasses when handling glass.

Consumables

Replacement consumables can be purchased under the following part numbers.

Part Number	Description
70-803-3009	DualSeal Ferrule (PKT. 5)
STT-16-100	Standard Sample Tubing 1/16" OD x 1.0mm ID x 10m long

For more information please contact one of our offices below.

Asia Pacific

6 Central Boulevard, Port Melbourne,
VIC 3207, Australia

Telephone: +61 3 9320 1111

Email: enquiries@geicp.com

www.geicp.com

Americas

31 Jonathan Bourne Drive, Unit 7 Pocasset,
MA 02559, USA

Telephone: 508 563 1800

Email: geusa@geicp.com

www.geicp.us

Europe

Friedenbachstrasse 9, 35781 Weilburg,
Germany

Telephone: +49 6471 3778517

Email: gegmbh@geicp.com

www.geicp.de

October 2026